

# Patented Maelstrom Spin Mixing Technology

TANBead Maelstrom product embodies this novel technology and delivers improved performance for applications in molecular diagnostics and life sciences. Maelstrom Series are FDA and CE approved, and the patents are granted in the Canada, China, EU, Korea, Japan, Taiwan, and USA.



#### **Fully Automated**

- Eliminate human errors , ensure quality consistency
- An open system for a broad range of magnetic bead-based reagent kits



#### Patented Whirl Stirring Mixing Technology

- Processing volume up to 1,500 μl
- Spin tips stir magnetic beads at speeds up to 3000 rpm



#### **Easy Operation**

- Intuitive user interface and easy menu navigation
- Data exchange via USB connected with laptop
- A portable device that can be used from lab to field



#### Time Saving

- 2,500 gauss magnetic rods efficiently collect magnetic beads
- Automation of complicated manual steps



# Operating procedures













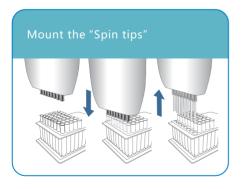



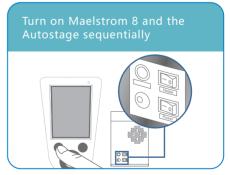







# Maelstrom<sup>TM</sup> 8 Autostage Series Features

## Performance

The coefficient of variation of nucleic acid extraction concentration is less than 5%



High consistency



No Cross-contamination to neighbor well Experimental results

### $1 \quad 2 \quad 3 \quad 4 \quad 5 \quad 6 \quad 7 \quad 8 \quad \mathsf{M} \quad 9 \quad 10 \ 11 \ 12 \ 13 \ 14 \ 15 \ 16$



# **Specification**


Model	Maelstrom™ 8
Run Time	25~60 min
Samples per run	1-8 samples
Weight (NW)	600 g
Dimensions (WxDxH)	11.2x6.3x32.7 cm
Power Supply	5 V , 2 A
Battery	3.7 V , 2,850 mAh
Charging Plug	micro-USB
Processing volume	50 μl - 1,500 μl
Magnetic Rod	Ø2.0mm, 2,500 gauss
Spin Speed	500~3,000 rpm
Display	2.4" TFT LCD

Model	M8-H Autostage
Weight (NW)	8.33 kg
Dimensions (WxDxH)	19x30.5x31 cm
Power Supply	19 V, 6.32 A 120 W
Temp control	1 set
Heating Block	Yes (1 pcs)
Heating	RT~110°C

## Patents

USA	US09616398B2	Korea	10-1696517
EU	EP2937136	China	CN104971638B
Canada	CA2862946	Taiwan	TWI526245B
Japan	JP6151735B2		

# Reagent kits

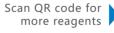





Sample	Target	Test	Kits Format	Ordering No.
Blood	Total DNA	96	Auto Plate	301128
Plant	Genomic DNA	96	Auto Plate	301134
Rice	Genomic DNA	48	Auto Tube	301137
cfDNA	cfDNA		Auto Tube	301155
FFPE	Total DNA	96	Auto Plate	301152
Viral	Viral OptiPure DNA / RNA		Auto Plate	301148
Tissue	e Total DNA		Auto Tube	301133
Fungi	Genomic DNA		Auto Plate	301150
Forensic	Forensic Total DNA		Auto Tube	301285
Gram Bacteria	Gram Bacteria Bacteria Total DNA		Auto Plate	301138


# Consumables

Images	Product Name	Format	Description	Ordering No.
	96 deep well plate	Auto Plate	<ul> <li>Processing volume 50 μl-1,600 μl</li> <li>Widely use for molecular diagnostics</li> </ul>	104070
	96-spin tips box	Auto Plate	• 96 pcs of Medium(pen-shape) spin tips in one box	303009
	6 tube B	Auto Tube	<ul> <li>Special Package for single or small number of tests</li> <li>Minimal consumable waste</li> <li>No reagent loss</li> </ul>	104143
	16-Base B	Auto Tube	Incorporate with 6 tube B for small number of tests	104026
	Spin tips	Auto Plate Auto Tube	A unique design to maximum mixing efficiency	104157

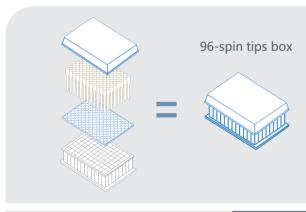


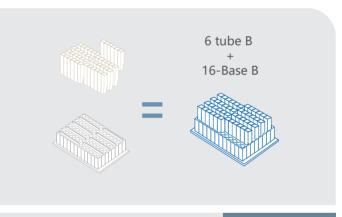



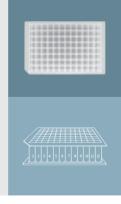



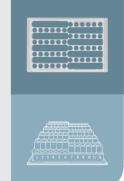



Taiwan Advanced Nanotech

www.tanbead.com